

NASA TECH BRIEF
Marshall Space Flight Center



NASA Tech Briefs announce new technology derived from the U.S. space program. They are issued to encourage commercial application. Tech Briefs are available on a subscription basis from the National Technical Information Service, Springfield, Virginia 22151. Requests for individual copies or questions relating to the Tech Brief program may be directed to the Technology Utilization Office, NASA, Code KT, Washington, D.C. 20546.


Mobile Automatic Metabolic Analyzer

The problem:

One of the goals in the design and fitting of prosthetics is to minimize the amount of energy that must be expended to use them. Energy consumption is normally measured while the device is being worn by a patient walking on a treadmill. In the case of above-the-knee amputees, however, treadmill exercise does not allow the forward motion necessary to operate the prosthetic.

The solution:

By adapting technology developed for NASA's Apollo program, a mobile automatic metabolic analyzer (MAMA) has been developed to measure and analyze breathed gases. The analyzer is housed in a mobile cart that overcomes the treadmill problem by automatically moving alongside a walking patient.

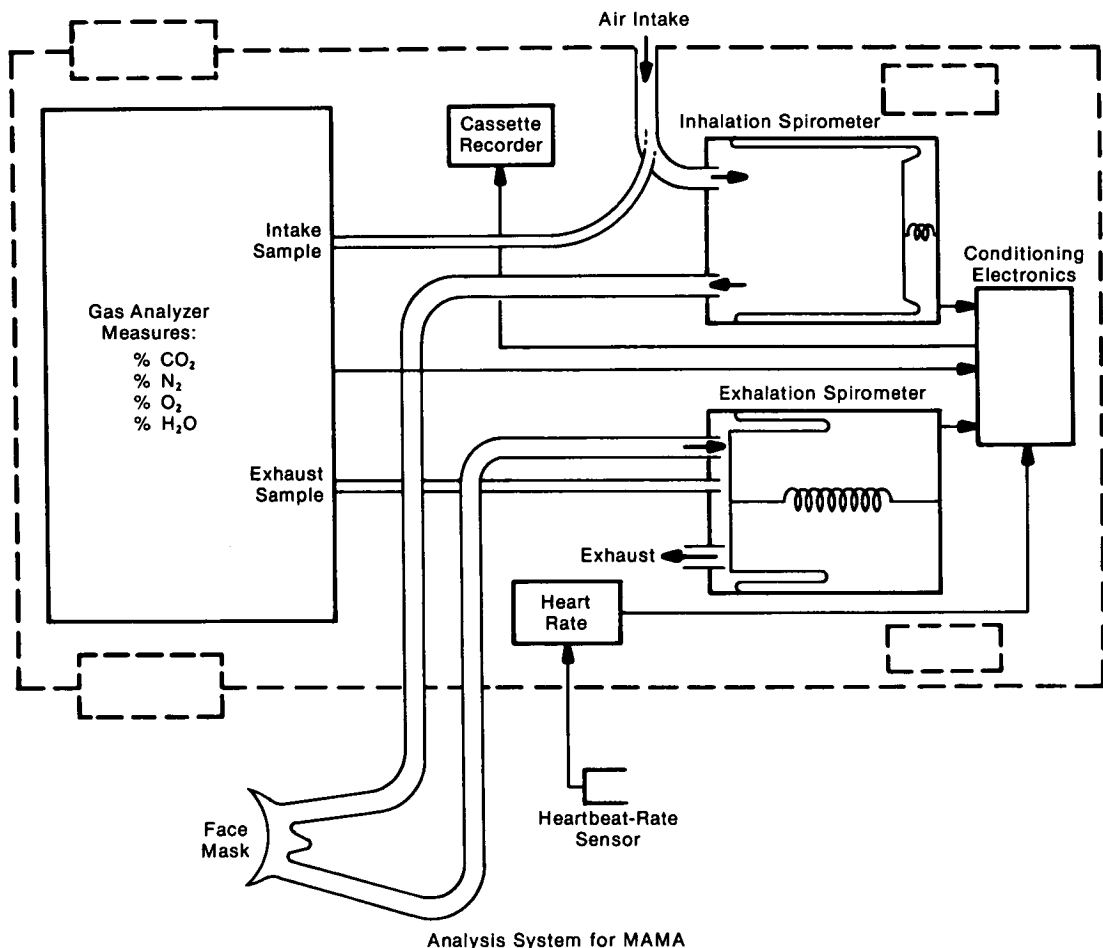

(continued overleaf)

How it's done:

The patient wears a face mask attached to two flexible pipes that are connected to spirometers in the cart. Inhaled air volume is measured as it is drawn from one spirometer, and exhaled air volume is measured as it is breathed into the second spirometer. A sensor, worn on the ear and connected to the cart, is used to monitor the patient's heartbeat rate. The medical equipment in the analyzer is shown in the illustration.

The heartbeat rate and the volumes of inhaled and exhaled air are each recorded on a tape as a function of time. Inhaled air is analyzed once prior to an experiment, and exhaled air is analyzed after each breath. These results are also stored on the tape.

(These data are later fed into a computer to determine the metabolic rate.) The analyzer also includes a handswitch on an extension cord, for use by the physician to mark events of interest on the tape, and to stop the cart.

The cart is powered by three battery-operated electric motors. It may be operated manually or automatically. In the automatic mode, the system optically senses a guidance tape laid out on the floor. The transistorized, automatic, speed and acceleration controls include safety features to accommodate possible patient problems. Some specifications of the mobile automatic metabolic analyzer are shown in the following table.

Acceleration	1 to 3 ft/s ²	(0.3 to 0.9 m/s ²)
Speed (Level Ground)	0 to 5.5 mi/h \pm 1%	(0 to 2.5 m/s)
Speed (10% Slope)	0 to 4 mi/h \pm 3%	(0 to 1.8 m/s)
Braking (Normal Stop)	4 ft/s ²	(0.9 m/s ²)
Braking (Emergency Stop)	5 ft/s ²	(1.5 m/s ²)
Turning Radius	3.8 ft	(1.16 m)
Operating Time	3 to 5 h	

Specifications for Mobile Automatic Metabolic Analyzer

Note:

Requests for further information may be directed to:

Technology Utilization Officer
Marshall Space Flight Center
Code AT01
Marshall Space Flight Center, Alabama 35812
Reference: B75-10077

Patent status:

Inquiries concerning rights for the commercial use of this invention should be addressed to:

Patent Counsel
Marshall Space Flight Center
Code CC01
Marshall Space Flight Center, Alabama 35812

Source: B. G. Bynum and J. R. Currie
Marshall Space Flight Center
(MFS-23143)